Expt. 22 Michelson interferometer

I. Purpose

- (1) Understand the theory, construction, adjustment and operation of the Michelson interferometer;
- (2) Learn the method of measurement of length using the Michelson interferometer.

II. Apparatus

Michelson interferometer, He-Ne laser, sodium lamp, cover glass, ground glass, observing panel with scales.

III. Principle

As shown in Figure 22 - 1 (a), a beam of parallel light travels from left to right to a beam-splitting lens. On passing through it is divided into two beams at the second surface of the beam-splitting lens (which is silver-gilt): One beam is reflected by the surface of the beam-splitting lens and then travels along one arm to the fixed plane mirror M_1 , and is reflected towards the beam-splitting lens and pass through it. The other beam travels along the other arm to the plane mirror M_2 , and is also reflected towards the beam-splitting lens and is reflected at the surface of it. The two beams of light interfere upon superposition and move towards the observation screen, forming beautiful interference fringes.


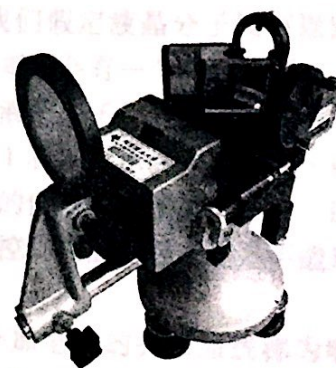
Figure 22 - 1 Michelson Interferometer

When the Michelson interferometer is used for measuring or calibrating one arm is always fixed. In this case M_2 is fixed and M_1 is moving. As M_1 is moved the interference pattern changes by

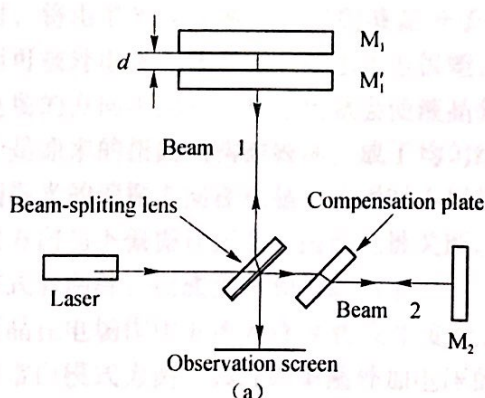N fringes and the displacement Δd of M_1 is given by

$$\Delta d = N \frac{\lambda}{2} \quad (22-1)$$

If the interference fringes are adjusted until they are clear and wide, then Δd can be readily measured from the movement of fringes.

IV. Experimental

1. Adjustment of the Michelson interferometer

(1) Turn on the He-Ne laser, set it approximately horizontal and make sure that the laser beam is parallel to the plane of the interferometer track. Adjust the three horizontal screws at the bottom of the interferometer to make sure that the laser beam reaches to the centre of M_1 and M_2 at the same time, after the beam has been split at the centre of the beam-splitting lens.

(2) Adjust the mirror adjustment screws on the back of M_1 and M_2 , to make sure that the central maximum light spots reflected from M_1 and M_2 arrive at the exit hole of the interferometer and superpose with each other.

(3) Once they superpose you can observe that the central maximum of the interference pattern appears black fringes. Keep on adjusting the screws of M_2 and do not move those of M_1 to make the fringes as wide as possible.

2. Observation of interference fringes of equal inclination

Adjust the height and position of the beam-expanding lens, so that the laser beam illuminates the entire beam-splitting lens P. Adjust in turn the "vertical adjusting screw" and the "horizontal adjusting screw" at the bottom of M_2 , so that you can observe the concentric annular fringes of equal inclination.

3. Observing the change of equal inclination fringes and analyzing the relationship between the change and the optical path difference of the two light beams

Adjusting the remote barrel wheel d , observe the fringes "disappearing" and "appearing" at the centre, and the change of width. Further more, you should determine the relationship between these phenomena and the optical path difference of the two light beams.

4. Measuring the wavelength of the He-Ne laser

The arm with the mirror M_2 is fixed while the other mirror M_1 is allowed to move along one direction. If you count the number of interference fringes N which pass a fixed point when the mirror M_1 moves from an original position d_1 to a final position d_N , you can find the displacement Δd , of M_1 . This is related to the wavelength λ by

$$\lambda = \frac{2\Delta d}{N} = \frac{2(d_N - d_1)}{N} \quad (22-2)$$

The procedure required is as follows:

Adjust the vertical and the horizontal adjusting screws under the mirror M_2 to move the annular interference fringes to the visual centre of the observation screen. Record the readings on the ruler scale, the barrel wheel D , and the remote barrel wheel d . Then turn the remote barrel wheel d

along one direction while 100 fringes appear or disappear from the central light spot. Record the new reading of the position of M_1 , d . Record 10 sets of data totalling 1000 fringes.

5. Observation of localized fringes

Localized fringes are obtained by adjusting the "mirror adjusting screws" on the back of M_2 and the vertical (or horizontal) remote screws and the remote barrel wheel, once the concentric annular fringes of equal inclination are observed. The interference fringes are observed to become wider and straighter. These are the well known localized interference fringes from the usual interference of light.

6. Observation of the interference of white light

As the white light comprises all colours with various wavelengths, besides the central nil-level bright fringe, the conditions for a bright fringe at a certain wavelength will coincide with the condition for a dark fringe at another wavelength. This makes it difficult to observe interference fringes for white light. To observe white light fringes the optical path difference between the two arms must be close to zero.

Determine the direction of movement of the barrel wheel D and remote barrel wheel d , when the fringes shrink and the central light spot disappears.

When the width of the fringes is about 7 – 8mm (about 1.5 divisions on the scales of the ground glass) and the fringes are comparatively straight then the optical path difference of the two arms is close to zero.

7. Measurement of the refractive index of a solid

The refractive index of a solid can be determined using the interference of white light. Here cover glass, which is a translucent and isotropic solid, will be used as the unknown sample. The schematic diagram of the optical path is shown in Figure 22 – 2.

In this experiment we use white light fringes. Firstly, adjust the apparatus and obtain white light fringes. Put the central dark fringe and coloured fringes on the visual centre of the observation screen. Record the position of M_1 . Then insert the cover glass with thickness D and refractive index n into the M_1 arm. Turn the barrel wheel D and the remote barrel wheel d anticlockwise, moving the mirror M_1 towards the beam-splitting lens P , decreasing the geometric distance between M_1 and P . Make sure the white light fringes appear at the visual centre again. Record the new position of M_1 . The displacement of M_1 , ΔL , is the optical path difference added by inserting the cover glass, and $\Delta L = (n - 1) D$, so we have

$$n = \frac{\Delta L}{D} + 1 \quad (22 - 3)$$

According to this formula, if the thickness D of the cover glass is known, the refractive index

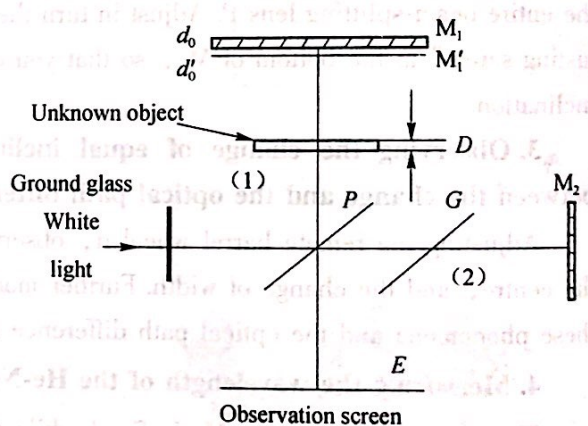
Figure 22 – 2 Determination of the Refractive Index of a Solid Using the Michelson Interferometer

of translucent solid can be worked out. Conversely, if the refractive index is known, the thickness can be determined.

V. Questions

1. How does the distance between the laser and the beam-splitting lens P affect the observation of the interference fringes of equal inclination? Does the position of the fringes of equal inclination change if the distance between the observation screen E and the beam-splitting lens P is changed? Does this affect the radius of the annular interference fringes? Try to explain all of your answers.

2. Why are white light fringes, which are difficult to obtain, used to measure the refractive index of the cover glass? Why do we not choose instead fringes of equal inclination, which are easily observed and measured?

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